




Polypropylene/Polyethylene Separator Depth Filters

Description

Virtually Inert to Chemical and Biological Agents

- ▶ Outstanding wet strength
- ▶ High chemical resistance
- ▶ Provides a channeling effect for liquid, increasing overall filter flow rate when used in series with submicron membrane filters
- ▶ Ensures stability when filtering concentrated acids, alkalis, or oxidizers

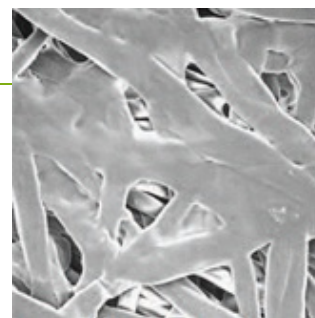

Application

- ▶ Good prefilter for heavily particulate-laden and viscous fluids
- ▶ Can be folded or handled the same as laboratory paper
- ▶ Can be used as separators between stacked discs in a single filter holder during serial filtration

Specifications

Filter Media

- ▶ Polypropylene/polyethylene

Pore Size

- ▶ 10 μm (nominal)

Typical Thickness

- ▶ 127 μm (5.0 mils)

Maximum Operating Temperature - Water

- ▶ 107 °C (225 °F)

Sterilization

- ▶ Provided non-sterile; autoclavable if desired

Related Products

- ▶ 13 mm Plastic Swinney Filter Holder
- ▶ 25 mm Filter Funnel, Stainless Steel
- ▶ 47 mm Filter Funnel, Stainless Steel
- ▶ 25 mm Easy Pressure Syringe Filter Holder, Delrin* Plastic
- ▶ 142 & 293 mm Stainless Steel Disc Filter Holders

Ordering Information

Part Number	Description	Pkg	Price	Qty
Polypropylene/Polyethylene Separator Depth Filters				
60343	10 μm , 25 mm	500/ pkg	NA	<input type="text" value="0"/>
60344	10 μm , 47 mm	100/ pkg	NA	<input type="text" value="0"/>

60346	10 µm, 127 mm	100/pkg	NA	<input type="text" value="0"/>
60341	10 µm, 142 mm	100/pkg	NA	<input type="text" value="0"/>
60345	10 µm, 265 mm	100/pkg	NA	<input type="text" value="0"/>
60342	10 µm, 293 mm	100/pkg	NA	<input type="text" value="0"/>

Contact Information

Pall Office(s)

World Headquarters

25 Harbor Park Drive
 Port Washington, NY 11050
 USA
 Phone: (516) 484-3600
 Alternate Phone: **1-800-289-7255**
 Fax: (516) 801-9754
[Driving Directions](#)
[Map](#)

Distributor(s)

VWR International

1310 Goshen Parkway
 West Chester, PA 19380
 USA
 Phone: 800.932.5000
 Phone: 610.431.1700
 Fax: 610.436.1760
 Web: <http://www.vwrsp.com/>

© 2014, Pall Corporation, Pall, and other names are trademarks of Pall Corporation.
 ® indicates a registered trademark in the USA.